



Advincula  
**DELINEATOR™**  
with Koh-Efficient® Technology

Advanced Uterine Manipulator  
Exceptional Strength and Performance  
In an Easy-to-Use Disposable

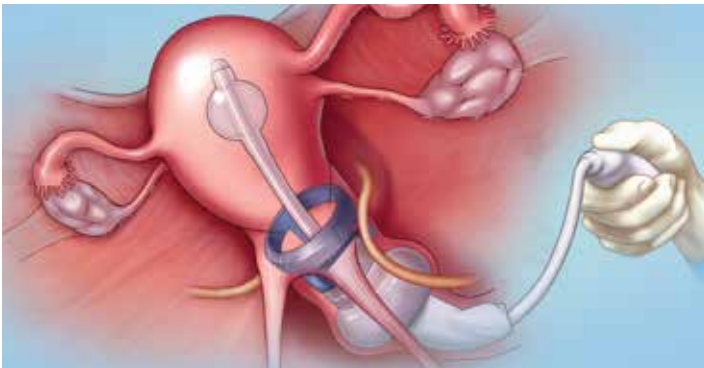


# Unprecedented Ease of Use, Safety And Strength for TLH, LSH and LAVH – In an Efficient Disposable Design

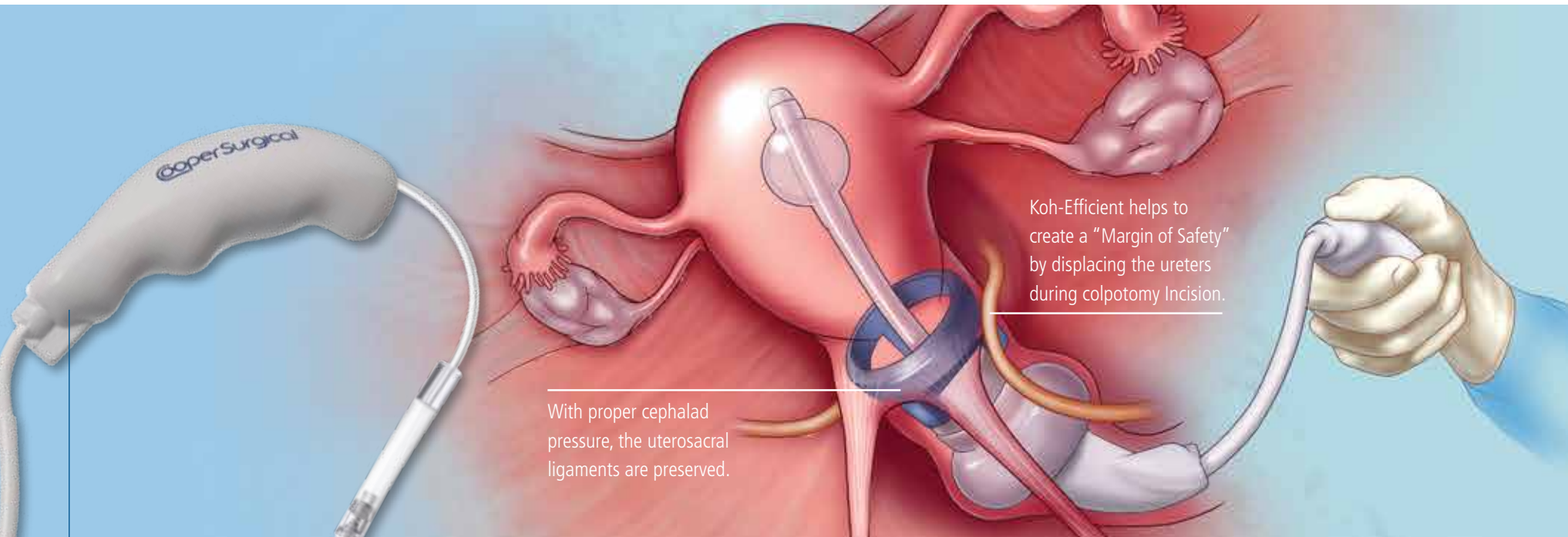

Koh-Efficient helps to create a "Margin of Safety" by displacing the ureters during colpotomy incision.

With proper cephalad pressure, the uterosacral ligaments are preserved.

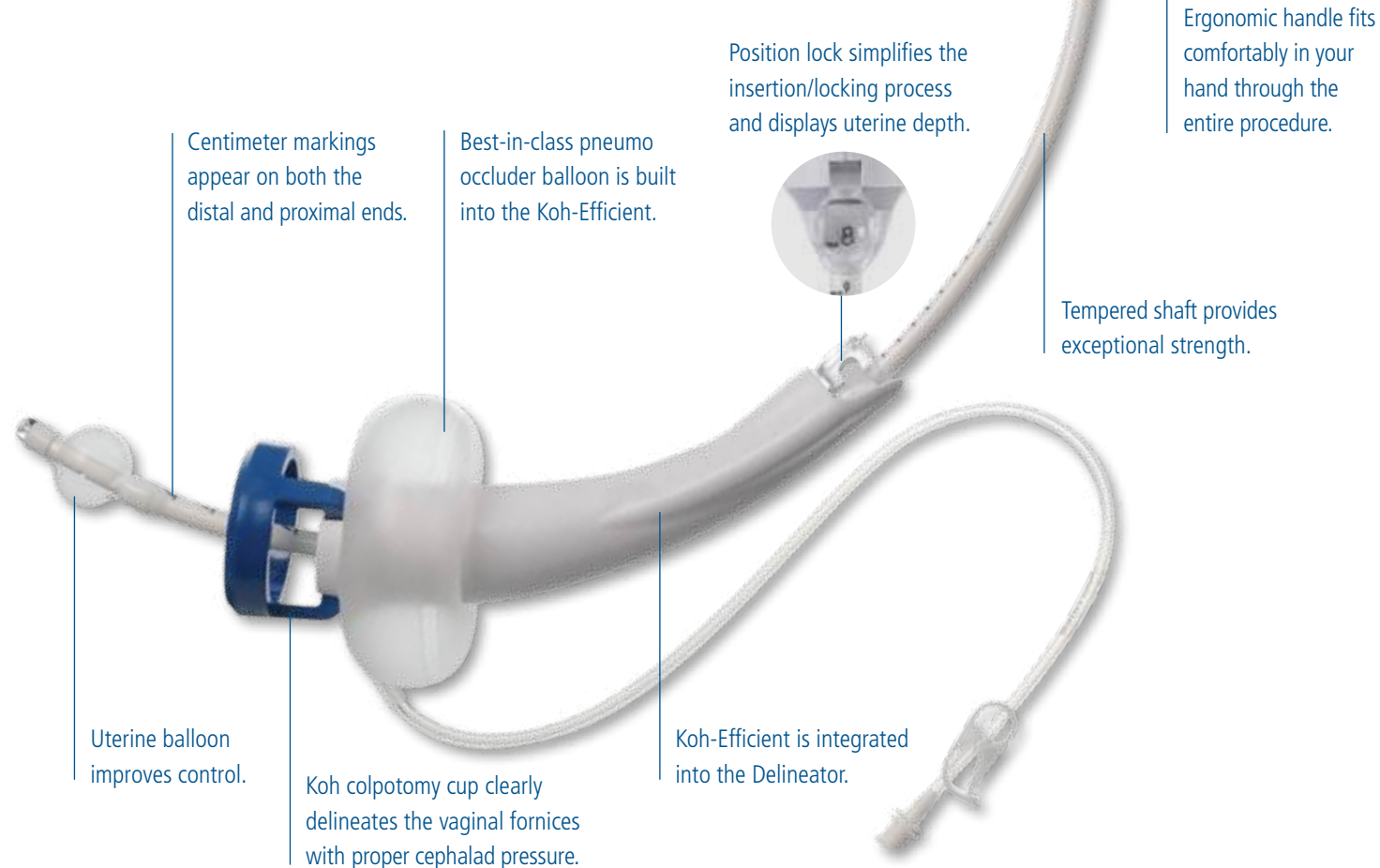

Centimeter markings appear on both the distal and proximal ends.

Best-in-class pneumo occluder balloon is built into the Koh-Efficient.

Position lock simplifies the insertion/locking process and displays uterine depth.



Tempered shaft provides exceptional strength.

Ergonomic handle fits comfortably in your hand through the entire procedure.

Uterine balloon improves control.

Koh colpotomy cup clearly delineates the vaginal fornices with proper cephalad pressure.

Koh-Efficient is integrated into the Delineator.

## Superior Design Simplifies Prep and Speed of Insertion

The Advincula Delineator™ is designed for ease of use, safety and maximum efficiency.


- Sound the uterus and use the Cervical Sizer to determine the appropriate Koh Cup size.
- Remove the fully-assembled device from the sterile packaging.
- Place manipulator through cervix and inflate uterine balloon.
- Slide the Koh-Efficient® to the sounded number marked on the handle.
- Press the position lock to secure the cup, and begin the procedure.



Integrated, unified design, simplifies setup and insertion.

## Koh-Efficient Helps Distance Critical Anatomy to Create a "Margin of Safety"

With the application of proper cephalad pressure, the Koh-Efficient helps to create a "Margin of Safety" by preserving the uterosacral ligaments and displacing the ureters during colpotomy incision.

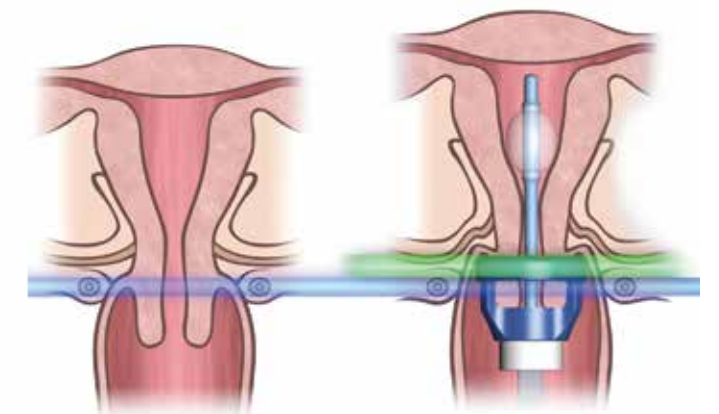

Available in  
Three Cup Options –  
To Accommodate  
Surgeon Preference  
and Energy Sources



Ultem™ plastic cup for use with  
electrosurgical devices.



Stainless steel cup for use with  
harmonic scalpel and laser devices.



Flexible, soft cup designed for easy insertion,  
and for use with electrosurgical devices.

# Advincula **DELINEATOR**™

with Koh-Efficient® Technology

## Ordering Information

### Advincula Delineator™ with Ultem Koh Efficient®

AD750-KE25	Advincula Delineator 2.5CM Ultem Plastic cup (3/Box)
AD750-KE30	Advincula Delineator 3.0CM Ultem Plastic cup (3/Box)
AD750-KE35	Advincula Delineator 3.5CM Ultem Plastic cup (3/Box)
AD750-KE40	Advincula Delineator 4.0CM Ultem Plastic cup (3/Box)

### Advincula Delineator™ with Stainless Steel Koh Efficient®

AD750-KES25	Advincula Delineator 2.5CM Stainless Steel cup (3/Box)
AD750-KES30	Advincula Delineator 3.0CM Stainless Steel cup (3/Box)
AD750-KES35	Advincula Delineator 3.5CM Stainless Steel cup (3/Box)
AD750-KES40	Advincula Delineator 4.0CM Stainless Steel cup (3/Box)

### Advincula Delineator™ with Flexible Soft Cup Koh Efficient®

AD750SC-KE25	Advincula Delineator 2.5CM Flexible Soft cup (3/Box)
AD750SC-KE30	Advincula Delineator 3.0CM Flexible Soft cup (3/Box)
AD750SC-KE35	Advincula Delineator 3.5CM Flexible Soft cup (3/Box)
AD750SC-KE40	Advincula Delineator 4.0CM Flexible Soft cup (3/Box)



For more information, call 800.243.2974 or 203.601.5200 or visit [www.coopersurgical.com](http://www.coopersurgical.com).  
800.243.2974 | 203.601.5200 | [www.coopersurgical.com](http://www.coopersurgical.com) | 95 Corporate Drive, Trumbull, CT 06611